


Quarter diopter IOL | Bi-Aspheric

Softec

HDY

High Definition Vision **Yellow**

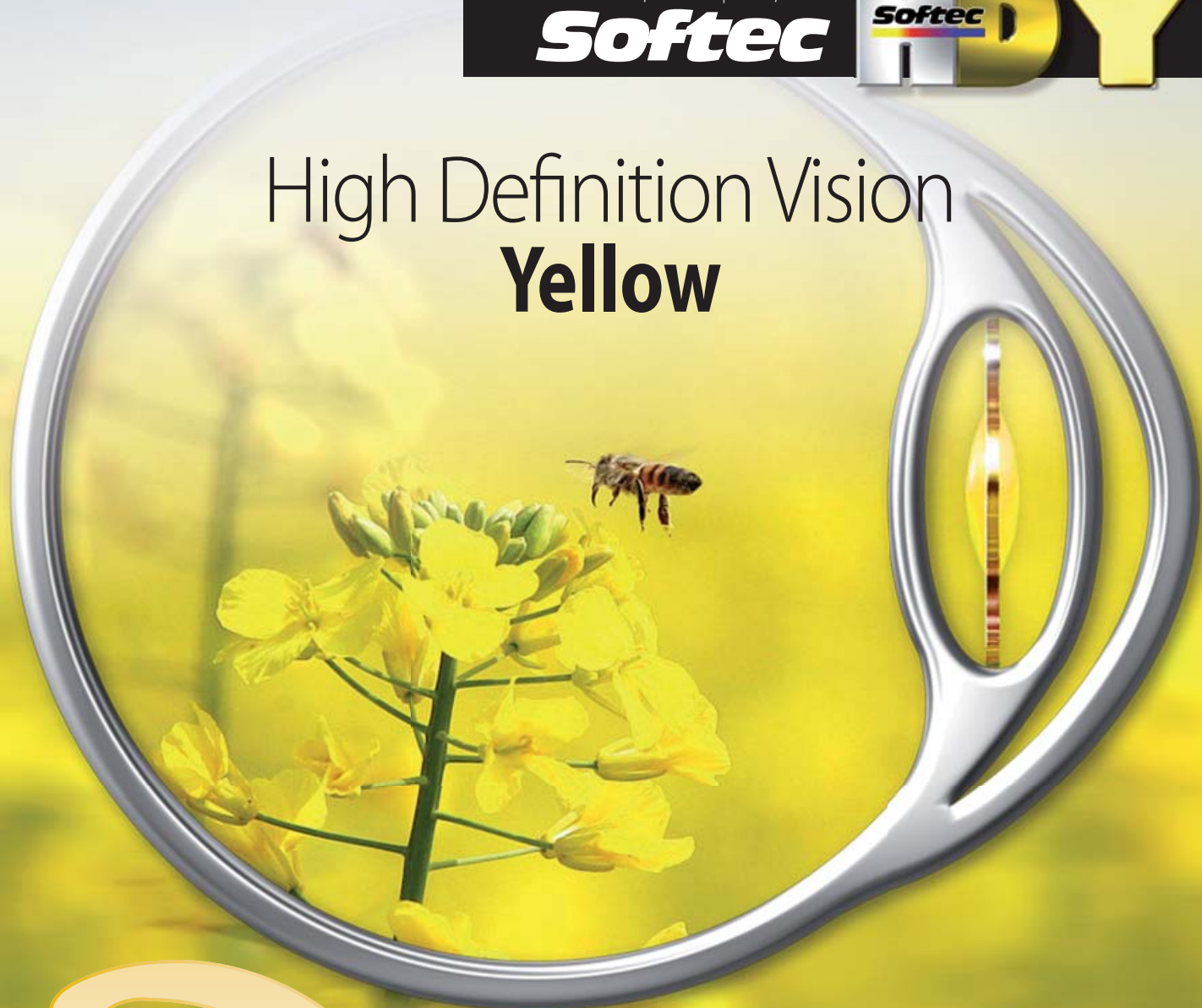

LENSTEGO™ LENSTEGO™ LENSTEGO™ LENSTEGO™ LENSTEGO™ LENSTEGO™

Single piece hydrophilic acrylic
Lenstec Softec HDY

- » Most accurate IOL with patented Bi-aspheric design with square edge technology
- » Greater precision due to quarter diopter increments
- » Tighter manufacturing diopter tolerance within ± 0.11 D
- » Proven quality for the most accurate cataract surgery

LENSTEC

more information on
www.lenstec.com

Technical Specifications

Optic Size	5.75 mm
Optic Type	Bi-aspheric
Length	12.00 mm
Haptic Style	Modified C
Angulation	0 Degrees
Positioning Holes	0
Construction	1 Piece
Optic Material	Acrylic (26% Water Content)
A/C Depth	5.10 mm
A-Constant(Contact Biometry)	118.00

A-Constant Optimized (Non-Contact Biometry)

Haigis	$a_0 = 0.92$ $a_1 = 0.40$ $a_2 = 0.10$
SRK/T	$A = 118.43$
Holladay1	$sf = 1.47$
Holladay2	5.22
Hoffer Q	$pACD = 5.22$

Diopter Steps

Whole	+5.00 to +36.00
Half	+10.50 to +29.50
Quarter	+15.00 to +25.00

Softec HDY - The Blue Blocking Lens

Single piece hydrophilic acrylic (26% Water Content)

The Softec HDY is fabricated from a material that has been proven stable, particularly with respect to long term exposure to sunlight, using an accelerated aging study that replicates the equivalent of 20 years of exposure for the average person (based on specific ISO standards for the photostability of IOLs). Using blue light absorbers in the fabrication of IOLs has been an ongoing discussion since the late 1980s, most researchers would agree that blocking certain wavelengths of visible light is beneficial to the eye... and a material proven to provide stable light transmission properties for over 20 years is the obvious choice.


Softec HD New Gold Standard - 0.11 D Tolerance - 3x More Precise!
Optical Prescription Selection and Tolerance Example: eye requires 24.25D Prescription to achieve optimal vision.

Industry Standard IOL			Softec HD IOL		
Lens Prescription	Allowed Tolerance		Lens Prescription	Allowed Tolerance	
24.00	±0.4	Max. Variance for a standard IOL in an eye that requires a 24.25D lens is 0.65D . (smaller number is better)	24.00	±0.11	Max. Variance for a Softec HD IOL in an eye that requires a 24.25D lens is 0.11D . (smaller number is better)
			24.25	±0.11	
24.50	±0.4		24.50	±0.11	
			24.75	±0.11	

Softec HD is the only IOL designed to address both Spherical Aberration and Defocus. Defocus is a more significant aberration than Spherical Aberration.

Bi-Aspheric Equal Conic Zero aberration IOL. Softec HD addresses the issue of spherical aberration inherent in conventional

monofocal spherical IOLs by adjusting the optic with a patented design on both the anterior and posterior surfaces.

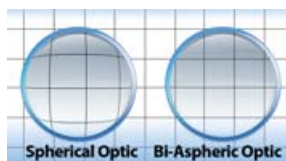

Studies have shown that Aspheric IOLs provide patients with significant optical benefits over traditional spherical surface IOLs.^{2,12,17}

Significant Outcomes. Mean refractive outcome was found to be closer to intended outcome, Depth of field was significantly improved, and Critical print size for fluent reading was smaller when compared to a standard monofocal IOL.⁷

Proven quality

Lenstec is one of eight companies in the world certified by FDA (Food and Drug Administration) for the sale of an intraocular lens in the U.S. market. All products have CE certificate, are approved by BSI (British Standards Institute) and are ISO quality system certified.

Stability of the biomaterial from which the intraocular Lenstec lenses are made, is long-term study proven and confirmed by millions of implanted lenses worldwide.

2. Thilos L, Hong X, Bradley A, Chang X. Statistical variation of aberration structure and image quality in a normal population of healthy eyes. J. Opt. Soc. Am A, Vol 19, No 21/Dec 2002 » 7. Craig J, Shah S, Wolffsohn J. Clinical evaluation of the Softec HD aberration-free aspheric intraocular lens. Submitted for publication. » 12. Sarver E. Theoretical optical performance of an equal conic intraocular lens and comparison to spherical and aspheric IOLs. AAO Presentation 2005 » 17. Nanavaty M, et al. Wavefront aberrations, depth of focus, and contrast sensitivity with aspheric and spherical intraocular lenses: fellow eye study. J Cataract Refract Surg. 2009; 35: 663 - 671

World Headquarters

Lenstec, Inc.
1765 Commerce Ave. N.
St. Petersburg, FL 33716
USA
Tel: 727-571-2272
Fax: 727-571-1792
Email: lenstec@lenstec.com

Lenstec UK

Lenstec (Barbados) Inc. Lenstec House
Unit 8, Mariner Court Calder Park
Wakefield WF4 3FL
UK
Tel: +44 (0)1924 382 678
Fax: +44 (0)1924 850 454
Email: lenstecuk@lenstec.com

Lenstec Barbados

Lenstec (Barbados) Inc.
Airport Commercial Centre
Pilgrim Road, Christ Church BB17092
Barbados
Tel: 246-420-6795
Fax: 246-420-6797
Email: lenstecbarbados@lenstec.com

Lenstec Regional Office

Trenčianska 47
821 09 Bratislava
Slovak Republic
cell: +421 905 798 760
e-mail: jkriska@lenstec.com